



ANALIZA BIAŁEK



100004129

## Laboratory of Analysis

**Client:** ep-company  
**Order number:** 100004129  
**Name:** EURO-PHARMACIES TEST 500  
**Number of samples:** 1 vial

**Batch number:** EP-4221C1  
**Recive date:** 2023-09-19  
**Start date:** 2023-10-02  
**End date:** 2023-10-03

## Sample identification

Substance	Method of taking a sample	Form
-----	Sample provided by customer	Transparent liquid in a glass vial


## Results

Substance	Sample name	Content	Purity
Testosterone Propionate	EURO-PHARMACIES TEST 500 500MG/ML	103,22mg/ml	-----
Testosterone Enanthate	EURO-PHARMACIES TEST 500 500MG/ML	190,48mg/ml	-----
Testosterone Cypionate	EURO-PHARMACIES TEST 500 500MG/ML	206,12mg/ml	-----

## Research comment

The test sample was applied in the amount of 3µl to the BioBasic 8 Thermo Scientific column. During the analysis with the method a linear gradient in the solutions of H2O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peaks of the sought substances appeared. Its concentration, when compared to the reference chromatogram shows the table above. The measurement error is + - 1mg. The actual state of the sample sent before the test is presented on the photo below

## Photo of the test sample

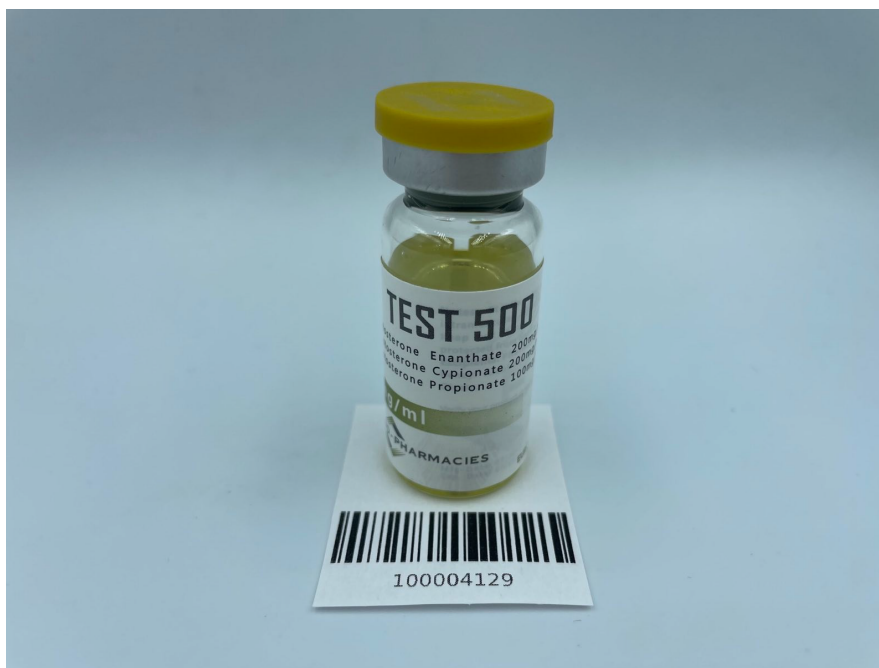

## Password

V37SJ032

**Analiza Białek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zajac




ANALIZA BIAŁEK



100004129

## Chromatogram

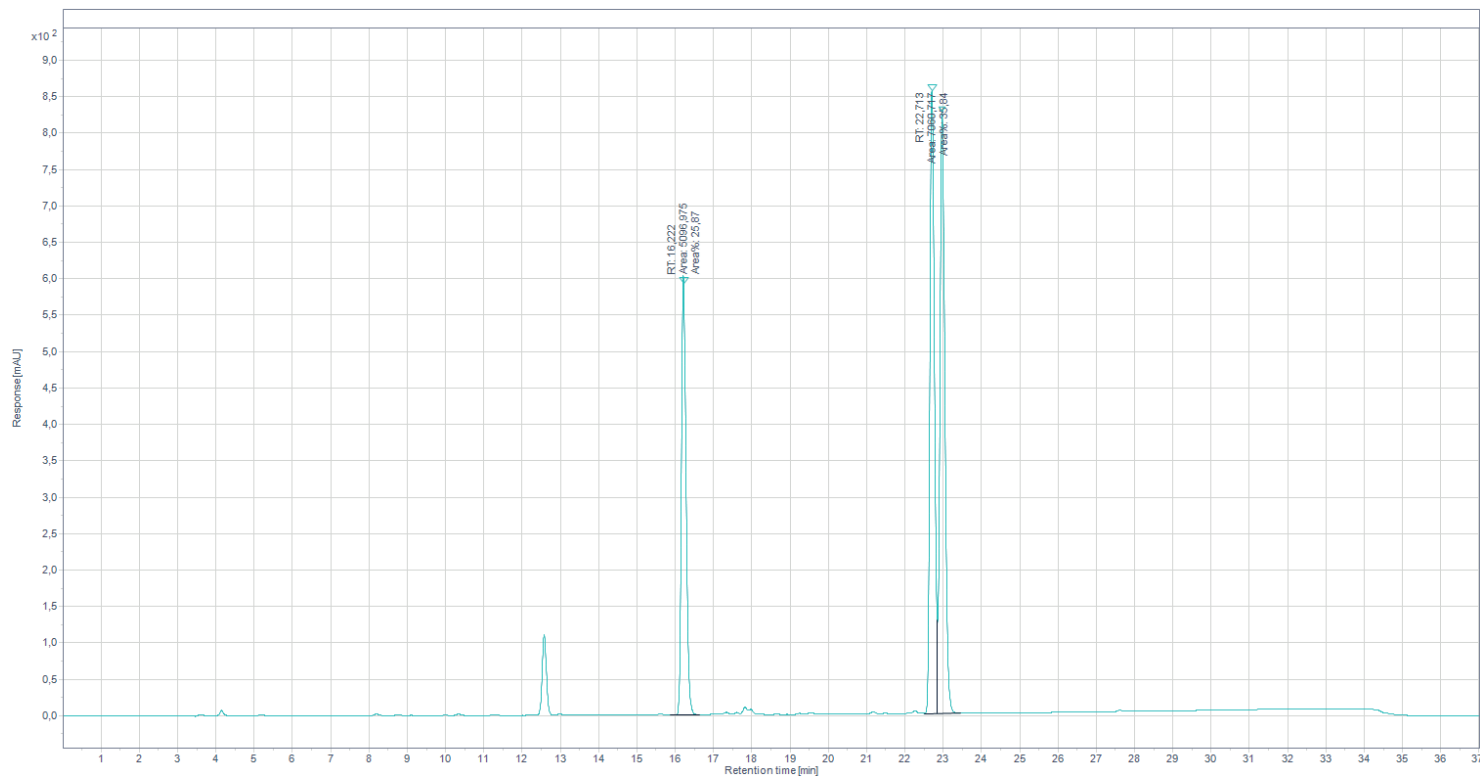

## Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

## Password

V37SJ032

**Analiza Białek sp. z o.o.**

**KRS:** 0000966168

**NIP:** 8992922970

**REGON:** 521751875

**Headquarters and Laboratory:**

ul. Duńska 9, 54-427 Wrocław

Budynek DELTA

**Lab technician:**

Analysis and research specialist

mgr Arkadiusz Zajęc