

ANALIZA BIAŁEK



100004127

## Laboratory of Analysis

**Client:** ep-company  
**Order number:** 100004127  
**Name:** EURO-PHARMACIES BLEND 450  
**Number of samples:** 1 vial

**Batch number:** EP-4221C1  
**Recive date:** 2023-09-19  
**Start date:** 2023-10-02  
**End date:** 2023-10-03

## Sample identification

Substance	Method of taking a sample	Form
Testosterone Enanthate Nandrolone Decanoate	Sample provided by customer	Transparent liquid in a glass vial


## Results

Substance	Sample name	Content	Purity
Testosterone Enanthate	EURO-PHARMACIES BLEND 450 450MG/ML	239,43mg/ml	-----
Nandrolone Decanoate	EURO-PHARMACIES BLEND 450 450MG/ML	203,72mg/ml	-----

## Research comment

The test sample was applied in the amount of 3 $\mu$ l to the BioBasic 8 Thermo Scientific column. During the analysis with the method a linear gradient in the solutions of H<sub>2</sub>O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peaks of the sought substances appeared. Its concentration, when compared to the reference chromatogram shows the table above. The measurement error is + - 1mg. The actual state of the sample sent before the test is presented on the photo below

## Photo of the test sample

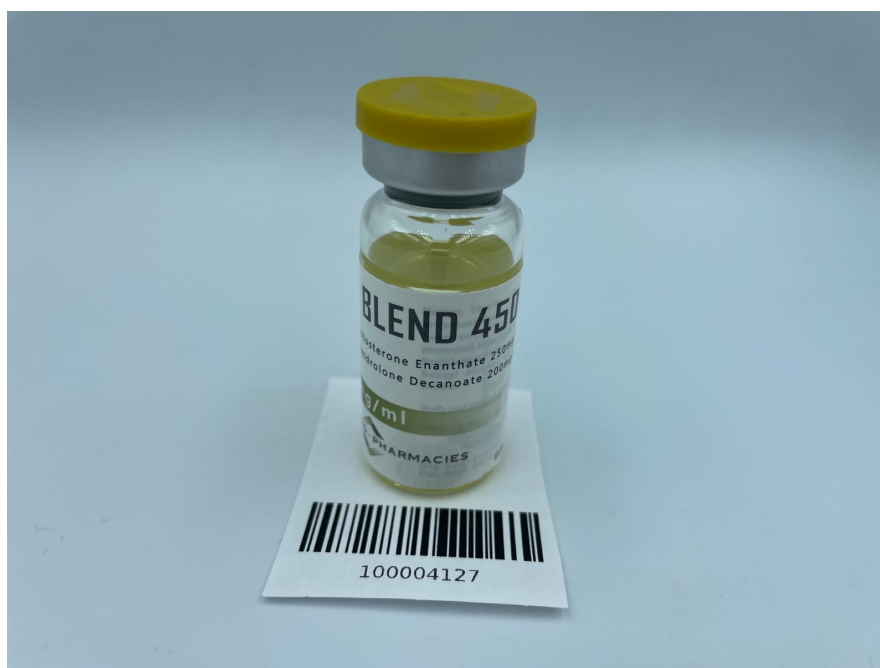

**Password**

N39BRD38


**Analiza Białek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zajęc



## Chromatogram

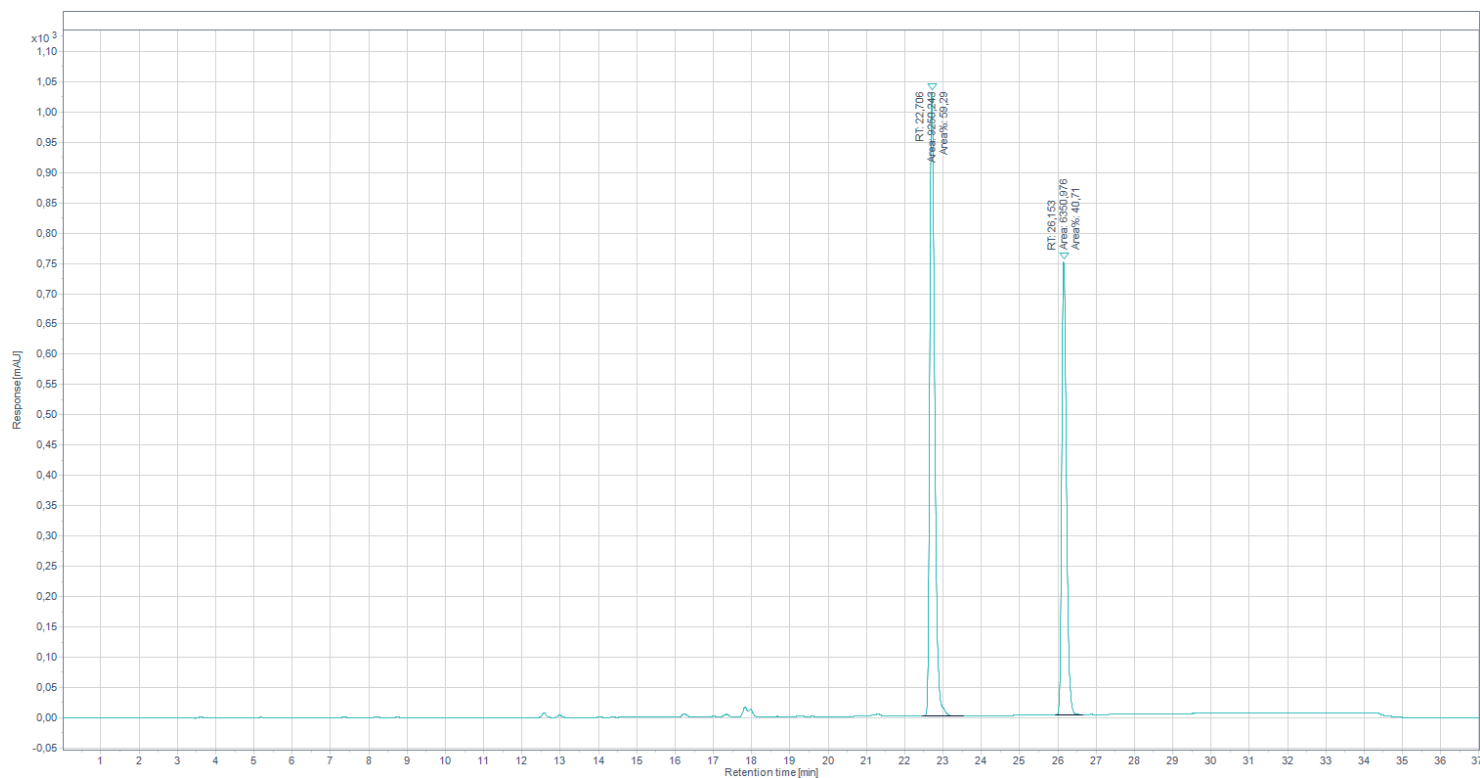

## Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

## Password

N39BRD38