



ANALIZA BIAŁEK



100005198

Laboratory of Analysis

Client: ep-company
Order number: 100005198
Name: EURO-PHARMACIES SUMMER CUTTING CYCLE
Number of samples: 1 tablet

Batch number: EP2200
Recive date: 2024-05-15
Start date: 2024-06-19
End date: 2024-06-26

Sample identification

Substance	Method of taking a sample	Form
-----	Sample provided by customer	White tablets in string bag

Results

Substance	Sample name	Content	Purity
Oxandrolone	EURO-PHARMACIES SUMMER CUTTING CYCLE 30MG+20MG	29,28mg/tab	-----
Stanozolol	EURO-PHARMACIES SUMMER CUTTING CYCLE 30MG+20MG	19,84mg/tab	-----

Research comment

The test sample was applied in the amount of 4 µl to the BioBasic 8 Thermo Scientific column. During the analysis with the method a linear gradient in the solutions of H2O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram shows the table above. The measurement error is + - 0,05mg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample



Password

2N7DHY43


Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajęc



Chromatogram

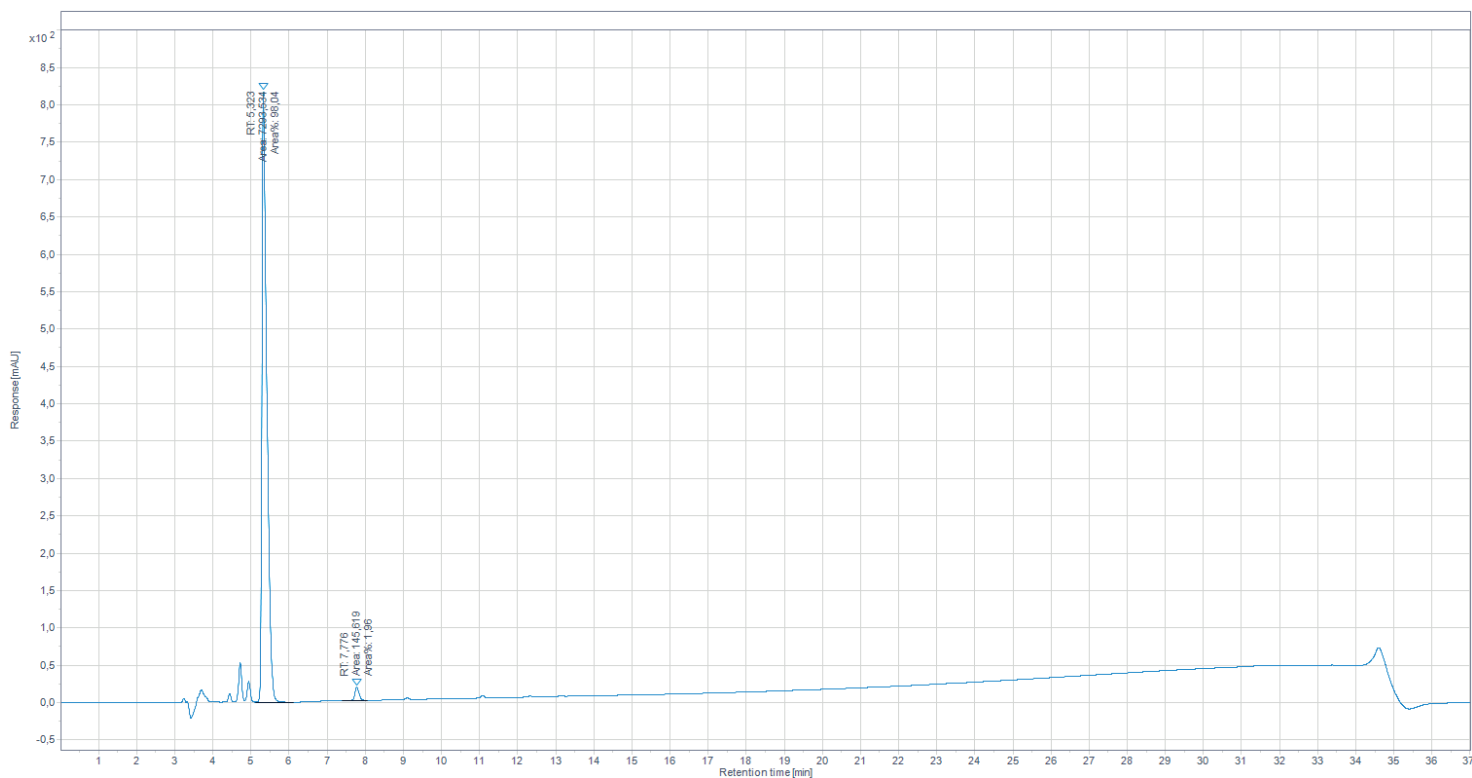

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

Password

2N7DHY43