



## **Laboratory of Analysis**

Client: ep-company Order number: 100005193 Name: EURO-PHARMACIES T4 Number of samples: 1 tablet **Batch number:** EP2500 **Recive date:** 2024-05-15 **Start date:** 2024-06-19 **End date:** 2024-06-26

## **Sample identification**

Substance	Method of taking a sample	Form
Levothyroxine	Sample provided by customer	White tablets in string bag

### **Results**

Substance	Sample name	Content	Purity
Levothyroxine	EURO-PHARMACIES T4 50MCG	46,13mcg/tab	

#### **Research comment**

The test sample was applied in the amount of  $10\,\mu$ l to the BioBasic 8 Thermo Scientific column. During the analysis with the method a linear gradient in the solutions of H20 and Acetonitrile with the addition of 0.1% trifuluoroacetic acid, a peak of the sought substance appeared Its concentration, when compared to the reference chromatogram, is 46,13mcg/tab. The measurement error is + - 0,5mcg. The actual state of the sample sent before the test is presented on the photo below

## Photo of the test sample




**Password** 

5B2M38KZ





# Chromatogram

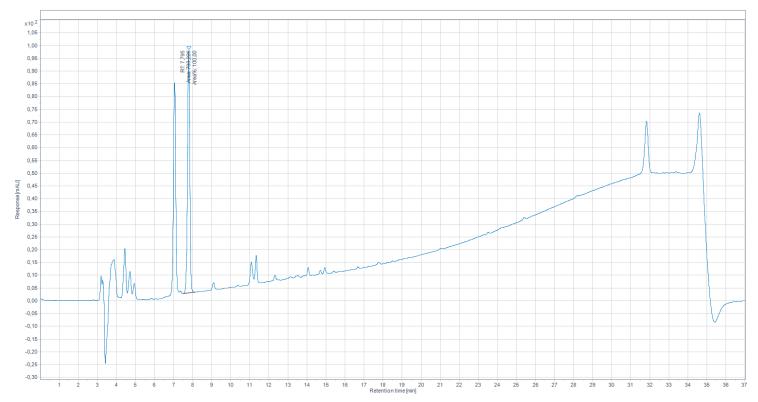

### **Method**

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

**Password** 

5B2M38KZ