



ANALIZA BIAŁEK



100009791

Laboratory of Analysis

Client: EP-COMPANY

Order number: 100009791

Name: EURO-PHARMACIES BOLDENONE UNDECYLENATE

Number of samples: 1 vial

Batch number: EP-Z7WDE3

Recive date: 2025-09-16

Start date: 2025-10-03

End date: 2025-10-06

Sample identification

Substance	Method of taking a sample	Form
Boldenone Undecylenate	Sample provided by customer	Transparent liquid in a glass vial


Results

Substance	Sample name	Content	Purity
Boldenone Undecylenate	EURO-PHARMACIES BOLDENONE UNDECYLENATE 250MG/ML	253,09mg/ml	-----

Research comment

The test sample was applied in the amount of 3µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H2O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 253,09mg/ml. Height the measurement error is + - 0,05mg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample

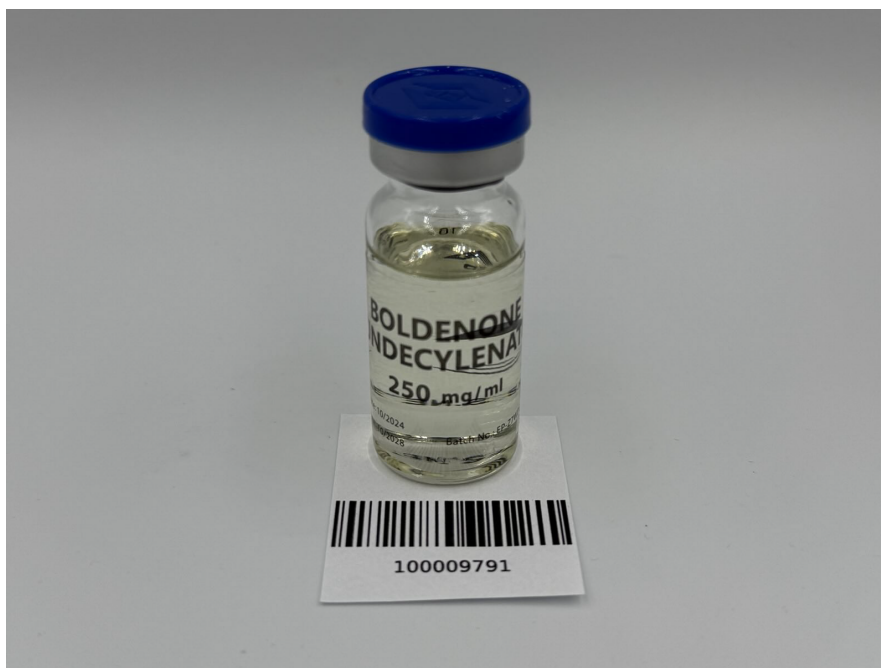

Password

3H7YR64V

Analiza Białek sp. z o.o.

KRS: 0000966168

NIP: 8992922970

REGON: 521751875

Headquarters and Laboratory:

ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:

Analysis and research specialist
mgr Arkadiusz Zajac




ANALIZA BIAŁEK



100009791

Chromatogram

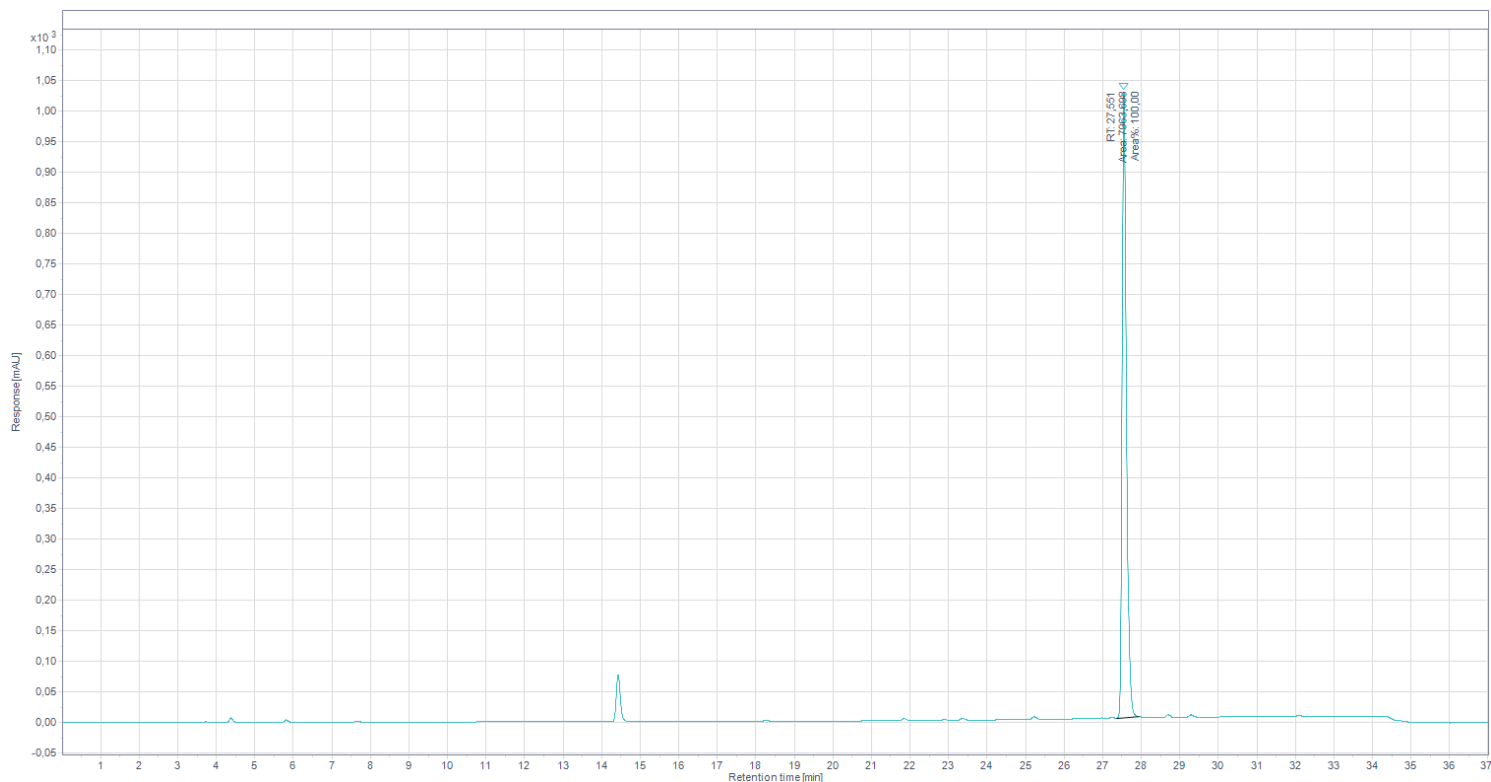

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

Password

3H7YR64V

Analiza Białek sp. z o.o.

KRS: 0000966168

NIP: 8992922970

REGON: 521751875

Headquarters and Laboratory:

ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:

Analysis and research specialist
mgr Arkadiusz Zajac