

ANALIZA BIAŁEK



100009813

Laboratory of Analysis

Client: EP-COMPANY
Order number: 100009813
Name: EURO-PHARMACIES EXEMESTANEX
Number of samples: 1 tablet

Batch number: EP1400
Recive date: 2025-09-16
Start date: 2025-10-06
End date: 2025-10-11

Sample identification

Substance	Method of taking a sample	Form
Exemestane	Sample provided by customer	White tablets in a string bag


Results

Substance	Sample name	Content	Purity
Exemestane	EURO-PHARMACIES EXEMESTANEX 20MG	21,09mg/tab	-----

Research comment

The test sample was applied in the amount of 3 μ l to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H₂O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 21,09mg/tab. Height the measurement error is + - 0,05 mg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample

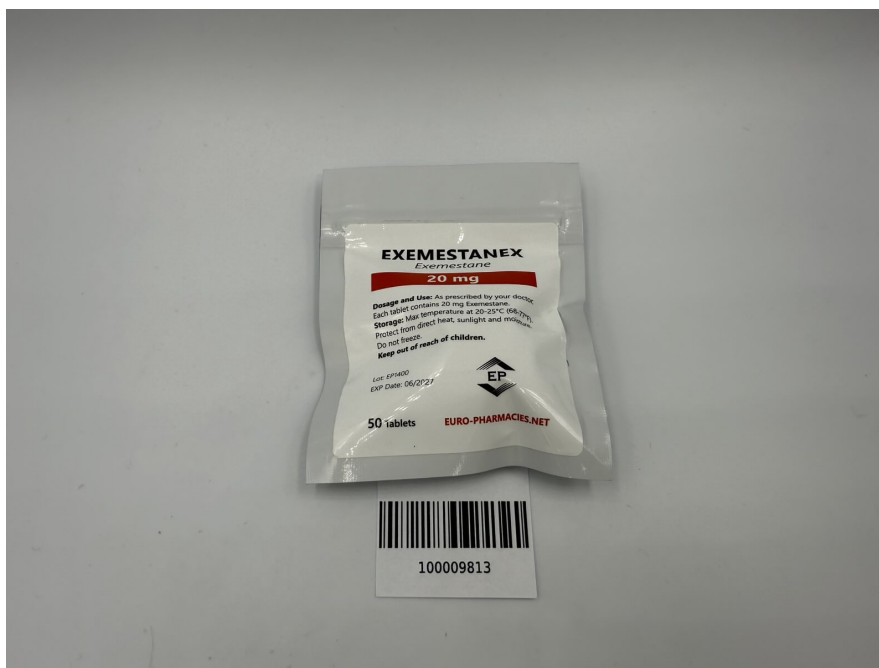


Password

G84LL9F3

Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajac

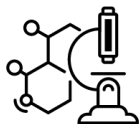


ANALIZA BIAŁEK



100009813

Chromatogram

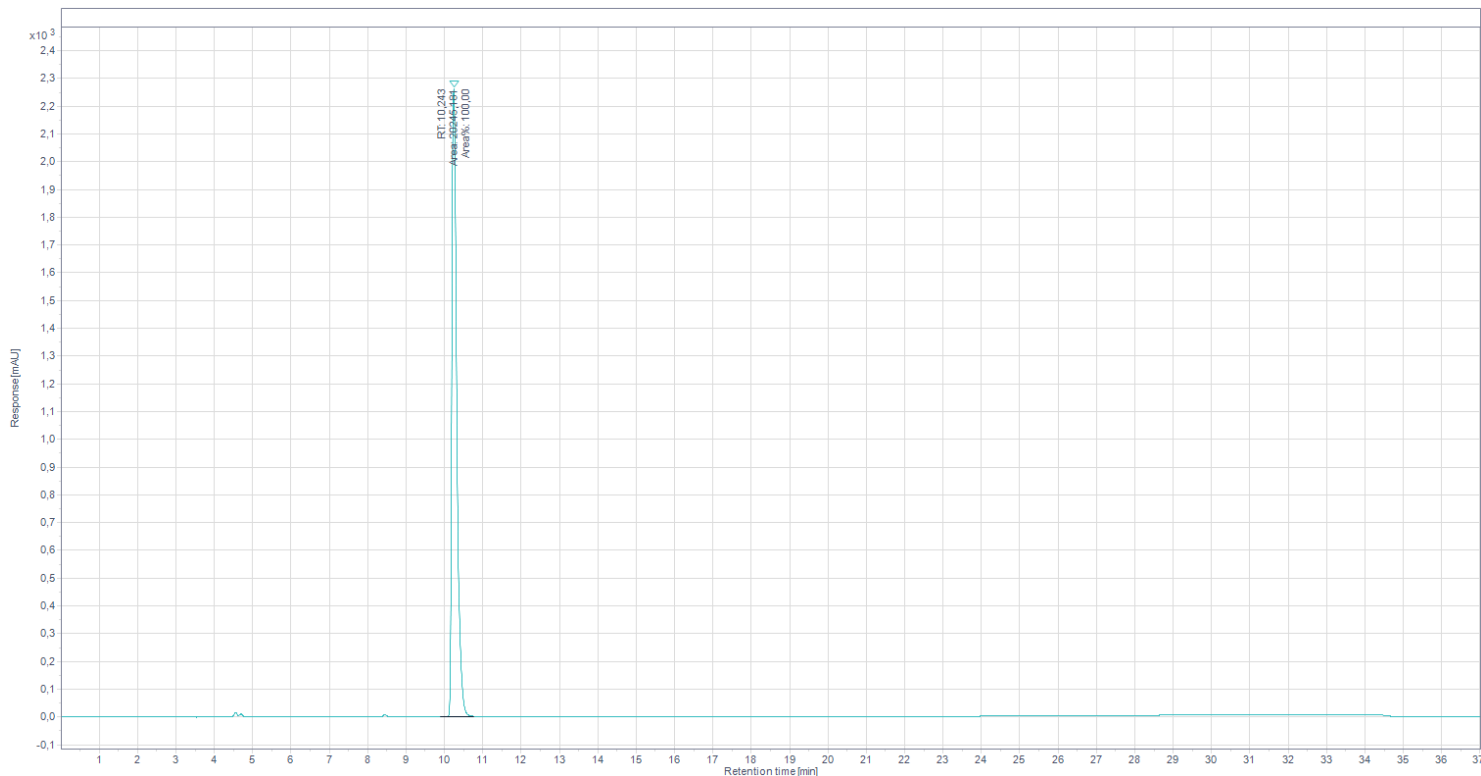

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

Password

G84LL9F3

Analiza Białek sp. z o.o.

KRS: 0000966168

NIP: 8992922970

REGON: 521751875

Headquarters and Laboratory:

ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:

Analysis and research specialist
mgr Arkadiusz Zajac