



ANALIZA BIAŁEK



100014678

Laboratory of Analysis

Client: EP-COMPANY
Order number: 100014678
Name: EURO-PHARMACIES.NET EXEMESTANEX
Number of samples: 1 tablet

Batch number: EP1400
Recive date: 2026-05-14
Start date: 2026-05-18
End date: 2026-05-20

Sample identification

Substance	Method of taking a sample	Form
Exemestane	Sample provided by customer	White tablets in a string bag

Results

Substance	Sample name	Content	Purity
Exemestane	EURO-PHARMACIES.NET EXEMESTANEX 20MG	19,65mg/tab	-----

Research comment

The test sample was applied in the amount of 2µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H2O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 19,65mg/tab. Height the measurement error is + - 0,05mg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample



Password

9K8U6J4V


Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajęc



Chromatogram

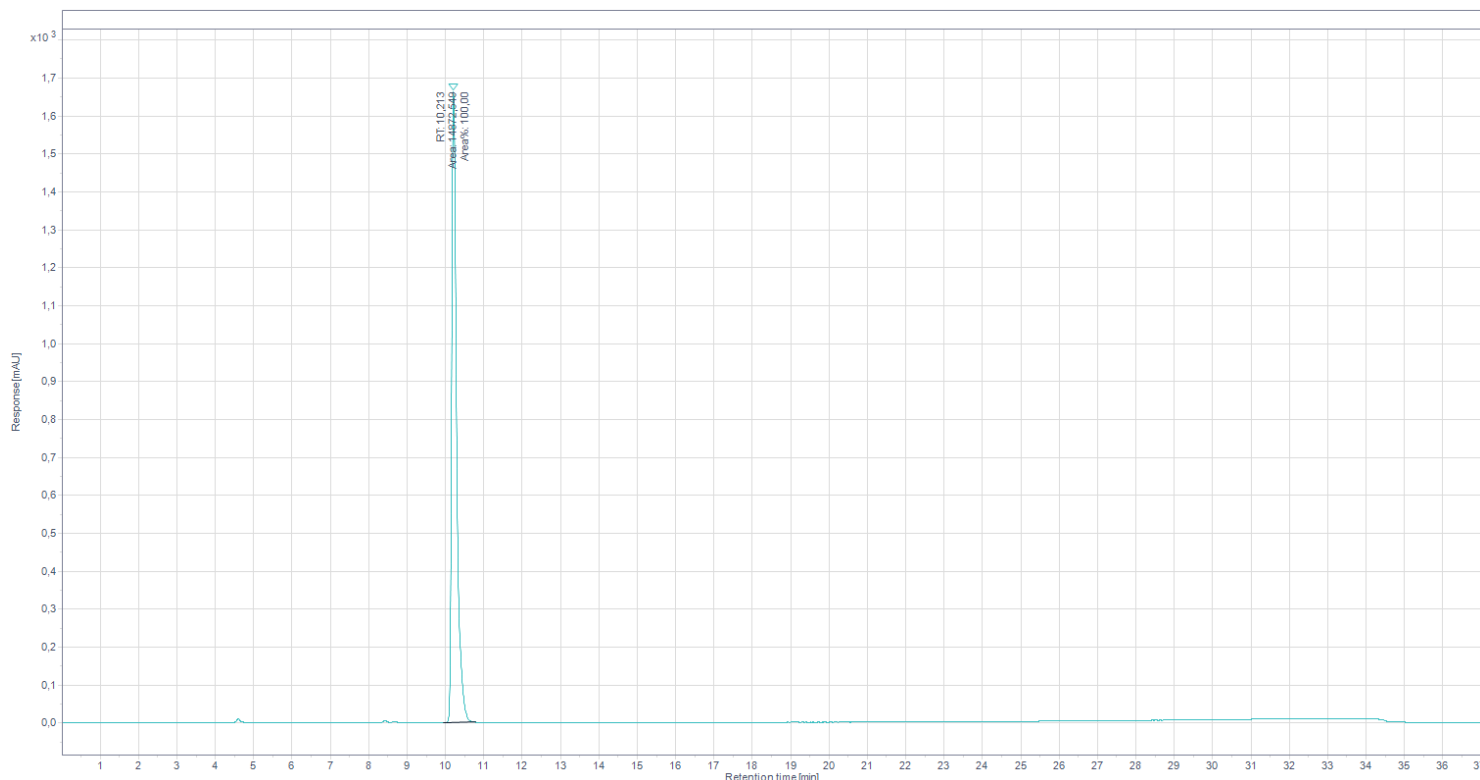

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	95%	5%
2min	95%	5%
12min	5%	95%
15min	5%	95%

Password

9K8U6J4V