

ANALIZA BIAŁEK



100015342

Laboratory of Analysis

Client: ep-company
Order number: 100015342
Name: EURO-PHARMACIES.NET CABASER
Number of samples: 1 tablet

Batch number: EP3100
Recive date: 2026-06-12
Start date: 2026-06-24
End date: 2026-06-26

Sample identification

Substance	Method of taking a sample	Form
Cabergoline	Sample provided by customer	White tablets in a string bag

Results

Substance	Sample name	Content	Purity
Cabergoline	EURO-PHARMACIES.NET CABASER 0.25MG	270,55mcg/tab	-----

Research comment

The test sample was applied in the amount of 10µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H₂O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 270,55mcg/tab. Height the measurement error is + - 0,5mcg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample




Password

8J74H3QN

Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajac

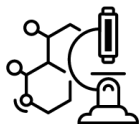


ANALIZA BIAŁEK



100015342

Chromatogram

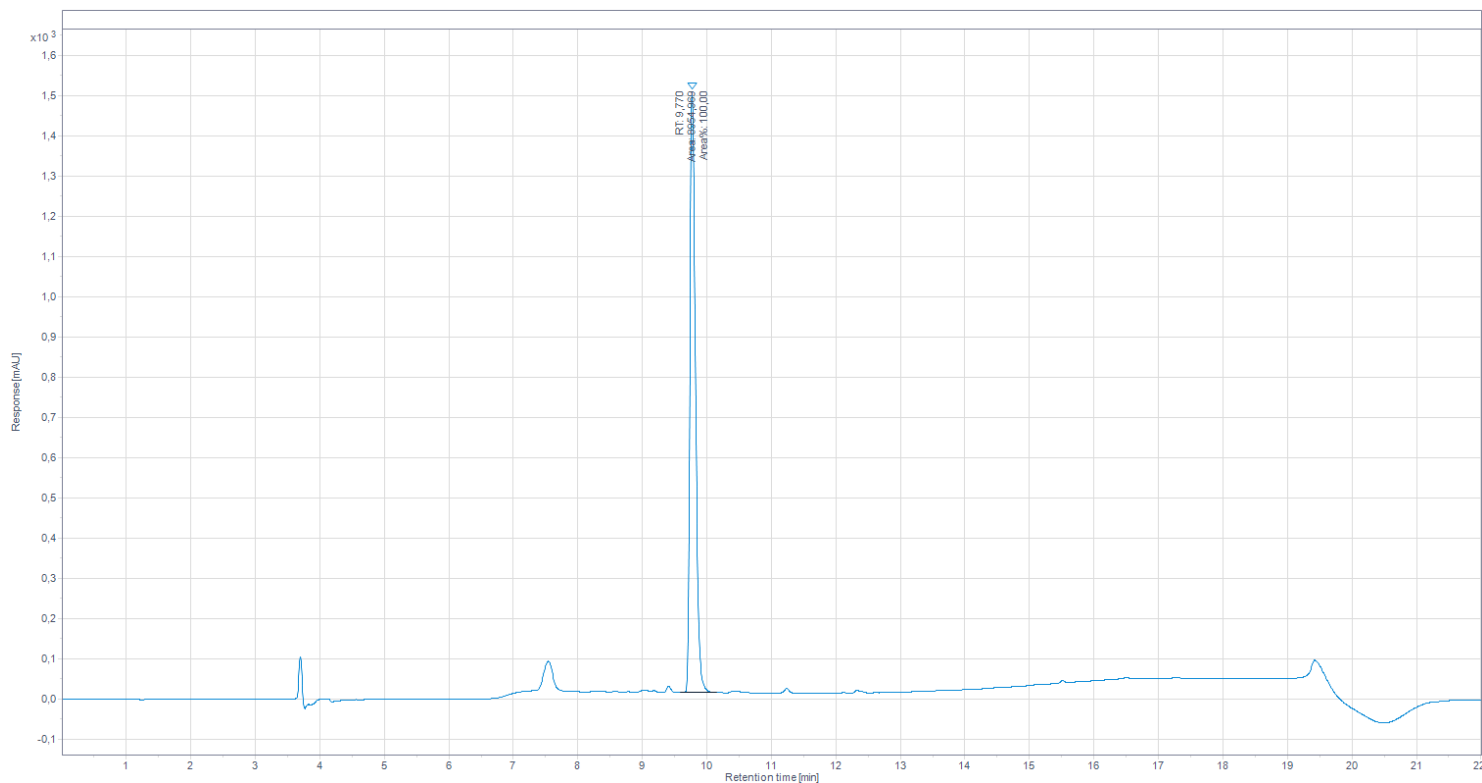

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	95%	5%
2min	95%	5%
12min	5%	95%
15min	5%	95%

Password

8J74H3QN

Analiza Białek sp. z o.o.

KRS: 0000966168

NIP: 8992922970

REGON: 521751875

Headquarters and Laboratory:

ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:

Analysis and research specialist
mgr Arkadiusz Zajac